


ITEM	QTY.	PART NO.	DESCRIPTION (MAIN)	DESCRIPTION (SUB)	MAT.
1	1	PA302300-	PARALLEL ACTION COMPRESSOR	AND DISTRACTOR SUB-ASSEMBLY	SEE CHART
2	1	PA302200-	HANDLE ASSEMBLY	DISTRACTOR	SEE CHART
3	1	PA202003	STANDARD BLADE SET	(5.5 - 6.35mm)	SEE CHART


REVISIONS			
REV	DESCRIPTION	DATE/BY	ECN NO.

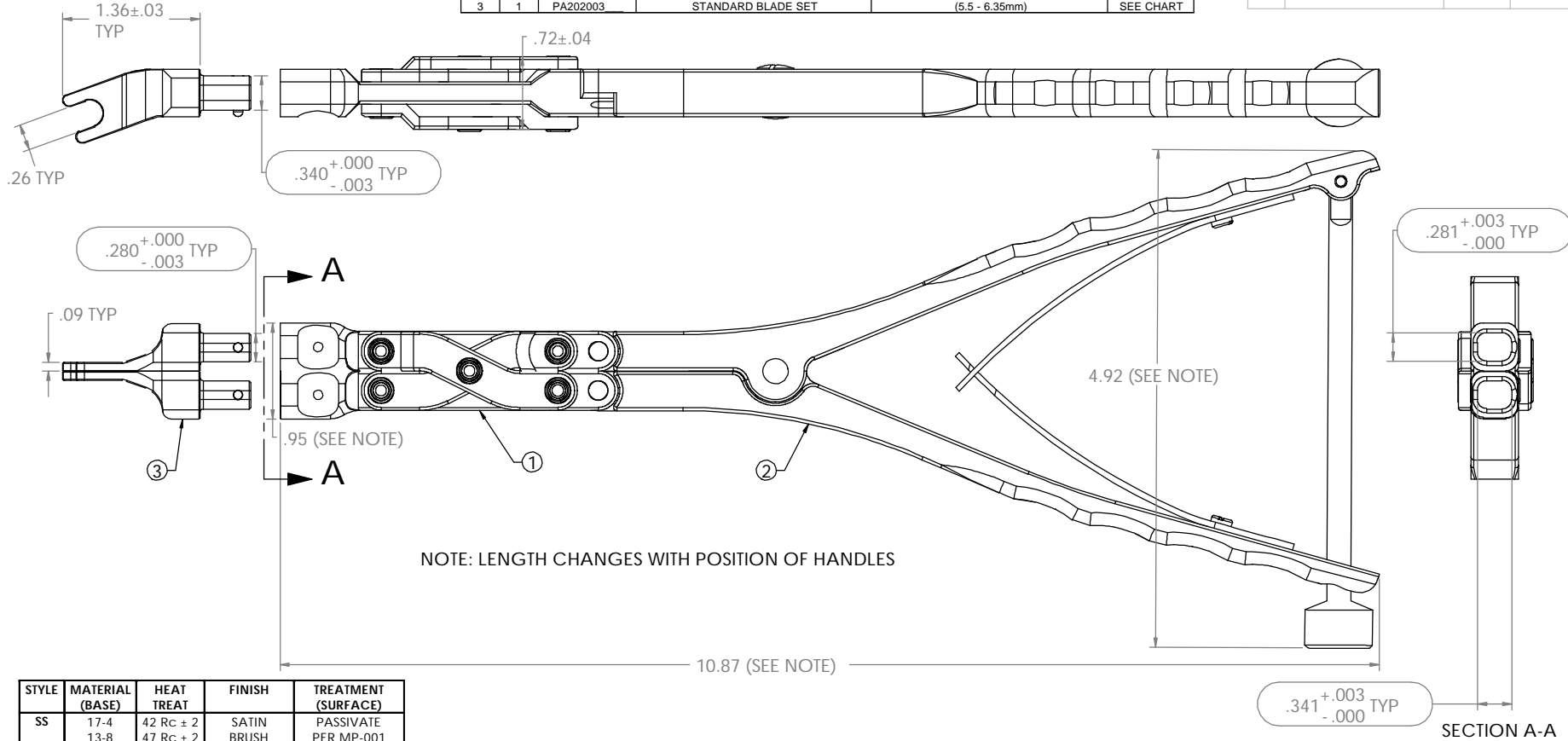

NOTE: LENGTH CHANGES WITH POSITION OF HANDLES

STYLE	MATERIAL (BASE)	HEAT TREAT	FINISH	TREATMENT (SURFACE)
SS	17-4 13-8	42 Rc ± 2 47 Rc ± 2	SATIN BRUSH	PASSIVATE PER MP-001
TI	6Al-4V	NONE	GLASS BEAD	PtFE
CR	6Al-4V	NONE	GLASS BEAD	CHROME PLATE
NT	6Al-4V	NONE	NONE	NONE

THIS DOCUMENT INCLUDING ALL DATA, SPECIFICATIONS, AND INFORMATION CONTAINED WITHIN IS CONFIDENTIAL AND PROPRIETARY TO BEERE PRECISION MEDICAL INSTRUMENTS, INC. BY RETENTION, THE RECIPIENT AGREES NOT TO REPRODUCE, DISCLOSE, OR MAKE USE OF ITS CONTENTS IN ANY WAY WITHOUT PRIOR WRITTEN PERMISSION OF AN OFFICER OF BEERE PRECISION MEDICAL INSTRUMENTS, INC.

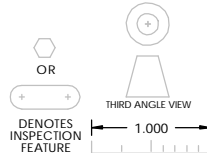

MATERIAL		UNSPECIFIED TOLERANCES		BEERE PRECISION MEDICAL INSTRUMENTS, INC.	
SEE CHART		.X +/- .020	ANGLE +/- .1°	2503 LINCOLNWOOD CT., RACINE, WI 53403 (262) 554-8505	
HEAT TREAT SEE CHART		.XX +/- .010	MACHINE FINISH 32	TITLE	
FINISH SEE CHART		.XXX +/- .005	CONCENTRICITY .005 T.I.R.	PARALLEL ACTION DISTRACTOR	
PA312203-		.XXXX +/- .0005	SQUARENESS .001 PER IN.	WITH INT. HANDLES & STD. BLADE SET	
DO NOT SCALE DRAWING		DRAWN JAR DATE 11/22/00		SIZE B	REV. 02
REMOVE ALL BURRS BREAK SHARP CORNERS		APPROVED DATE		SUPERSEDED NO.	DRAWING NUMBER PA312203-
PART TO BE FREE OF DIRT, OIL & CONTAMINANTS		SCALE 1:1		ILLUSTRATION SHEET 1 OF 1	